


# Brilliant images at every treatment unit

FLAT SCREEN AND ACCESSORIES FROM DÜRR DENTAL



TFT monitor installed at the treatment unit

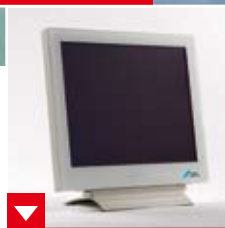

TFT monitor with  
foot stand



Easy adaptation of the  
handpiece holder and  
base unit



Different installation  
solutions

## Dürr Dental TFT monitor

17" TFT monitor with aluminum housing and protective glass with an optical-interference anti-reflective coating disinfected with wipe disinfectants, integrated foot stand and DVI-connection for digital data transfer

Ideally suited for: VistaCam PC and Video solution as well as VistaScan and VistaRay applications

- ▶ VistaCam handpiece holder can be mounted on the right of the monitor housing
- ▶ VistaCam base unit can be mounted at the rear of the TFT monitor

# TFT monitors for medical use

## PRACTICAL INSTALLATION SOLUTIONS FOR DÜRR DENTAL TFT MONITOR

### Optional accessories

- ▶ Different range of arm systems for wall fixing or fixing to the bar holding the lamp
- ▶ Cable set for connecting the VistaCam base unit to the TFT monitor
- ▶ Shared fastening rings for fixing to the bar holding lamp available in Ø 45, 50, 55 and 60 mm



TFT monitor

### Installation solutions for TFT monitors



Single-jointed arm



Double-jointed arm



Fastening ring



Wall mounting plate

Technical data	17" TFT monitor
Brightness	270 cd/m <sup>2</sup>
Contrast ratio	450 : 1
Resolution	SXGA 1280 x 1024
Connections	DVI-D, VGA, Video, S-Video
Installation	Mounting facility for VistaCam base unit in the rear of the TFT monitor and handpiece holder on the right of the monitor housing. TFT monitor is mounted via VESA Standard 100 Adapter on a single or double-jointed arm
Accessories	inclusive of medical grade power supply and 7 m extension cable for medical grade power supply.
Dimensions	408 x 398 x 55 mm (WxHxD)
Weight	6.8 kg
Color	RAL 9002 (as jointed arms for Dürr Dental TFT monitor)
Test mark	CE, TÜV GS complies with ISO 13406/2 in combination with VistaCam as per EC directive 93/42 EEC (Medical products law)

Dürr Dental GmbH & Co. KG  
 Höpfigheimer Strasse 17  
 74321 Bietigheim-Bissingen

www.duerr.de  
 info@duerr.de

