




Veraview iX Intraoral X-ray Compact Design with Exceptional Image Quality

- Easy to position with compact head design
- Small control panel minimizes installation space
- Patient protection with constant voltage control
- LED buttons with intuitive graphics
- Easy to switch between film and digital timer tables

Thinking ahead. Focused on life.



Equipped with LED guide for focal point positioning

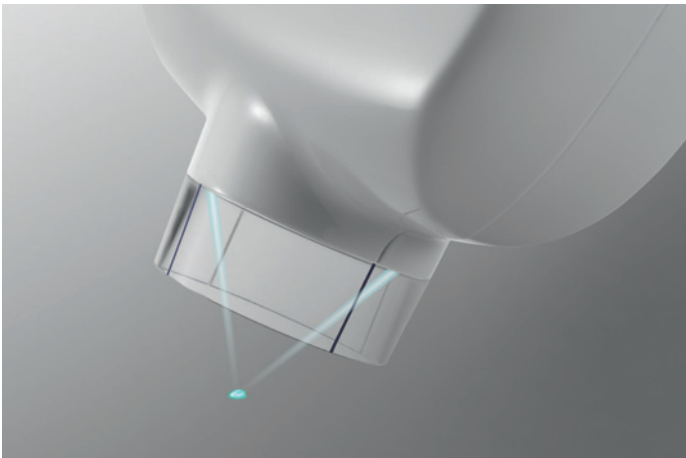

Quickly confirm position with the transparent cone and LED guide



- **Ease of use positioning features**
 - Positioning improved with the lightweight and compact design from pivot to cone.
 - Transparent cone to visually confirm position.
 - Easy to understand LED guide
 - Slender, lightweight arm for smooth, effortless movement
- **Small control panel minimizes installation space**
 - Small and compact, can be installed virtually anywhere.
- **Patient protection with constant voltage control**
 - Constant voltage control minimizes patient soft x-ray absorption.
 - Tube voltage options include 60 or 70 kV.
- **Easy to switch between film and digital timer tables**
 - Settings automatically saved in memory for both film and digital timer tables.



- **LED buttons with intuitive graphics**
 - See a visual representation of the current selection in the control panel for quick confirmation.

Specifications

Model	Veraview iX
Voltage	100-120 / 220-240 VAC
Tube Potential	60 / 70 kV
Tube Current	7mA
Exposure Time	10 ms to 2 seconds
Weight	approx. 23 kg

Developed and Manufactured by:
J. Morita Mfg. Corporation
 680 Higashihama Minami-cho, Fushimi-ku, Kyoto,
 612-8533 Japan
 T +81 (0)75. 611 2141, F +81 (0)75. 622 4595
www.jmorita-mfg.com

Distributed by
J. Morita Corporation
 33-18, 3-Chome, Tarumi-cho Suita City, Osaka, 564-8650 Japan
 T +81 (0)6. 6380 1521, F +81 (0)6. 6380 0585
www.morita.com/asia

J. Morita Europe GmbH
 Justus-von-Liebig-Str. 27a, 63128 Dietzenbach, Germany
 T +49 (0)6074. 836 0, F +49 (0)6074. 836 299
www.morita.com/europe

J. Morita Corporation Australia & New Zealand
 Suite 2.05, 247 Coward Street, Mascot, NSW 2020, Australia
 T +61 (0)2. 9667 3555, F +61 (0)2. 9667 3577
www.morita.com/anz

J. Morita Middle East
 4 Tag Al Aoasaa, Saba Pacha 21311, Alexandria, Egypt
 T +203 (0)58. 222 94, F +203 (0)58. 222 96
www.morita.com/middleeast